



Pure-Flo[®] & EnviZion[®] High Purity Diaphragm Valves

Product Selection Guide



ITT

ENGINEERED FOR LIFE

Section C

Integrated Block Valves

The Pure-Flo IBV - Integrated Block Valve product line is available in standard Pure-Flo and EnviZion body types. The valves are constantly changing to meet the rigorous demands of the Biopharmaceutical processing industry. An extensive array of innovative integrated block valves specifically designed to achieve the utmost in process efficiency. Pure-Flo developed the first integrated block body diaphragm valves over 30 years ago. We have a history of listening to customers and industry needs to develop valve solutions for the toughest applications.

Biopharm processes are complex and sensitive to system and environmental factors. Drug purity and process yield is greatly affected by system design. Integrated block valves can play a substantial role in developing a robust high yield process. Many valve solutions can produce acceptable results, but Pure-Flo integrated block valve technology can make a marginal process better and a good process great.


Integrated block technology is a cost effective means of reducing total cost of ownership. By optimizing drainability, hold-up volume, deadlegs and cleanability, block technology can decrease cleaning cycle times and increase process efficiency. Combining multiple valves into a single valve body can substantially reduce total installation and validation costs. Efficient designs pay for themselves over and over again.

Utilizing powerful 3D modeling software we can create almost any valve configuration imaginable. Working hand in hand our engineers will develop the valve configuration that fit your needs to a "T".

Table of Contents

Block Valve Benefits	C2-3
6D Rule vs. ASME BPD L/D	C4
P&ID Cross Reference	C5-6
Zero Static Use Points	C7-11
Divert & Sterile Access Valves	C12-17
Vessel Valves	C18-20

Integrated Block Valves in a Typical Bioreactor Process

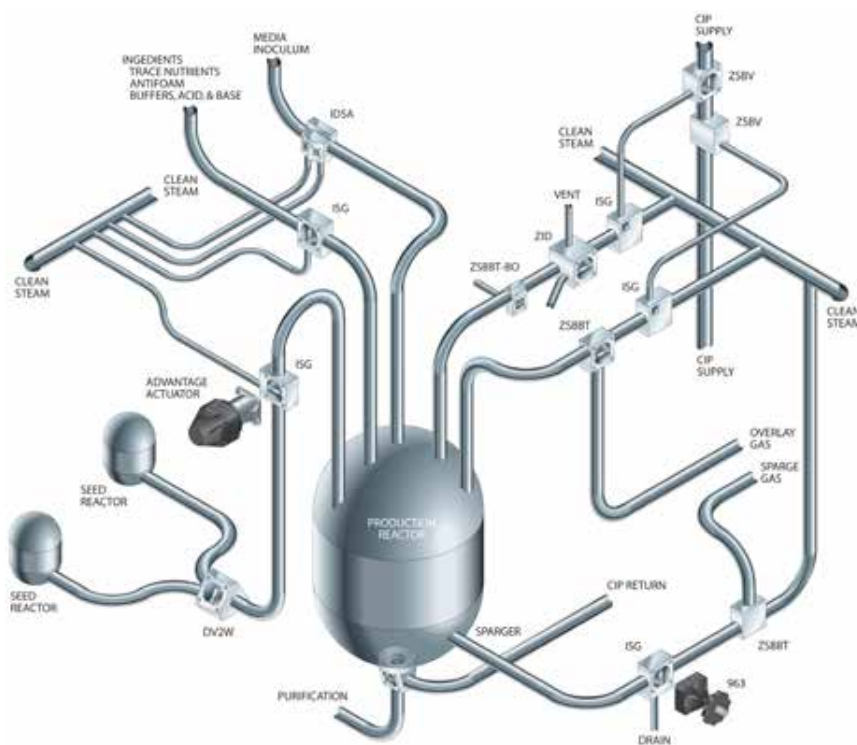

Block Valve Benefits

Drainability and Hold Up Volume

ITT Pure-Flo hygienic weir style diaphragm valves have become the most important control element of process piping systems utilized in the Pharmaceutical and Bioprocessing industries. The weir style diaphragm valve has become the standard due to its unique ability to provide maximum drainability and minimized product hold up volumes.

Integral Block technology further improves drainability and minimizes hold up volumes by reducing the process pipe volume between control elements.



Utilizing the unique characteristics of the weir style diaphragm valve, valve manufacturers have helped develop many process fabrications that have reduced product contact surfaces, reduced hold up volumes, and minimized piping dead-legs. The theory is that as contact surfaces are minimized, and hold-up volume are reduced in a process piping system, product yields and product purity will be improved.

Deadlegs

The FDA Guidelines for High Purity Water Systems has "defined dead legs as not having an unused portion greater in length than six diameters of the unused pipe, measured from the axis of the pipe in use. It should be pointed out that this was developed for hot (75-80° C) circulating systems. With colder systems (68-75° C) any drop or unused portion of any length of piping should be eliminated if possible, or have special sanitizing procedures".

In the not so distant past, typical process fabrications were produced by welding of standard forged valve bodies in configurations designed specifically for certain applications and orientations. This fabrication has served the industry well, but has limitations. In many instances the dead leg between the two can fall outside FDA expectations.

4-Way Divert Valve vs. Conventional Divert Valve Assembly

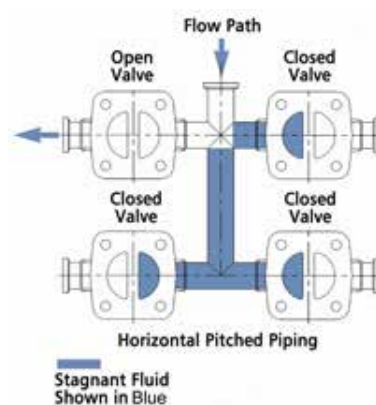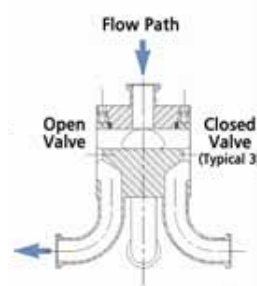

In case where process piping falls outside of FDA expectations, as noted from the High Purity Water Guide reference above, the owner of the system is expected to have special sanitizing procedures. These special sanitizing procedures can be costly in production time and processing cost and should be avoided whenever possible.

Block Valve Benefits

Current Good Manufacturing Practice (cGMP)

The cGMP regulation is a total quality concept applicable to processes and associated operations that assure the desired quality product. cGMP compliance, like quality, is fundamental and must be designed and built in from the earliest stages of a drug production project.

Drug manufacturers are required to maintain current Good Manufacturing Practices. This means that manufacturers must stay current with:

- New Technology
- New Methodology
- New Thinking
- New Requirements
- New Trends

Block Valves: Total Cost of Ownership

Total cost of ownership for a process system can not be calculated by material costs alone. Installation and ongoing operational costs should be taken into account when making any component purchasing decision. In many cases the cost of integrated block valves are greatly offset by reductions in installation costs, space requirements and improvements in operational efficiency.



One of the most critical factors in the production of drugs is the ability to clean and validate the drug production process. cGMPs require that processing equipment be designed to be cleaned and sterilized to minimize the potential for contamination, assuring the purity of the end drug product.

Hygienic weir style diaphragm valves have become the most important control element of process piping systems utilized in the Pharmaceutical and Bioprocessing industries, due to their unique ability to provide drainability with minimized product entrapment areas. Integrated block valve designs take these characteristics to an even higher level.

Integrated Block Valves can improve production efficiencies by:

- Minimizing internal valve volume
- Minimizing hold up
- Minimizing dead-legs
- Reducing CIP cycle times
- Increasing product purity
- Reducing qualification and validation efforts

Integrated Block Valves also reduce:

- Installation time and cost
- Expensive field welds
- Process piping footprint



ENGINEERED FOR LIFE

ITT Engineered Valves
33 Centerville Road
Lancaster, PA 17603, USA
Tel: +1 (717) 509-2200

**Cam-Line, Cam-Tite, Dia-Flo,
EnviZion, Pure-Flo, Skotch**

ITT Engineered Valves
1110 Bankhead Avenue
Amory, MS 38821, USA
Tel: +1 (662) 256-7185

Fabri-Valve

ITT Industries Ltd.
Weycroft Avenue,
Millwey Rise Industrial Estate
Axminster, EX13 5HU, United Kingdom
Tel: +44 1297-639100

EnviZion, Pure-Flo